




Engineered for Excellence :

Western Blot Imaging Meets Semiconductor-Grade Technology

Chip size

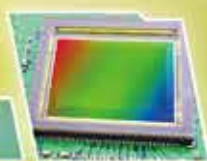

Western blot



Signal-Direct Imager

Chemiluminescence Western Blot Imaging System





Leveraging Taiwan's world-class semiconductor industry and its rapid advancement in production inspection equipment, SisBio Limited partners with local optoelectronics firms to transfer high-precision optical inspection technology directly into a superior Western Blot solution. This, combined with our 30 years of protein research expertise, delivers a cost-effective, high-resolution chemiluminescence system




Signal-Direct Imager

Redefining Western Blot Imaging Standards


-  Powerful simplicity
-  See more, 100x more sensitive
-  Measure everything, 100x wider range
-  A Footprint as Small as a Laptop
-  Get results in seconds

See more, 100x more sensitive.

This innovative design eliminates many bottlenecks of traditional Cooled CCD camera imaging

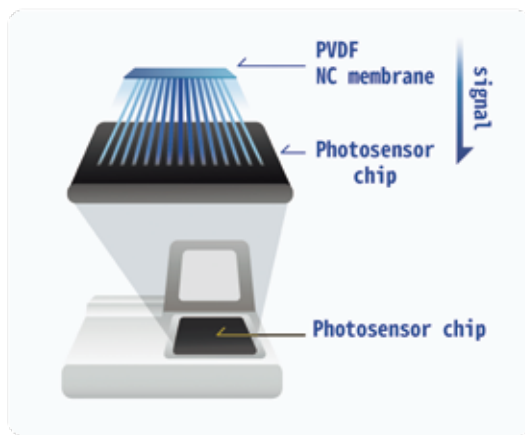

Signal-Direct Imager

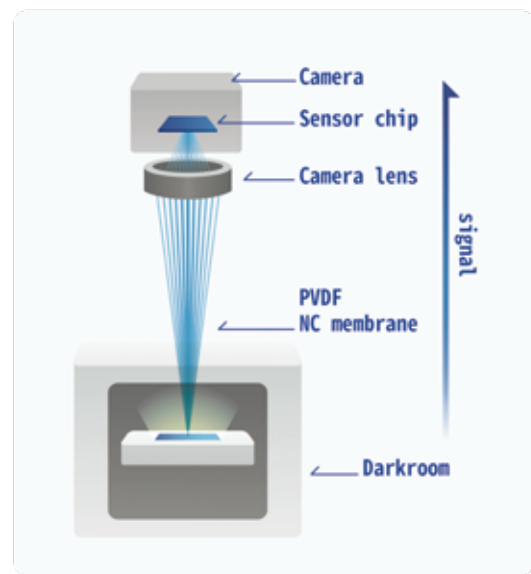


Cooled CCD camera imager

The sensitivity of a western blot imager is primarily determined by its ability to collect as many photons as possible from the blot and convert them into a strong, quantifiable signal while minimizing noise.

As shown in the image on the left, the new imaging technology removes any lens for imaging, thus eliminates image distortion caused by the lens.

The western film is directly attached to the photosensor chip for imaging, all signals are captured without any light loss, which greatly improves the sensitivity and speed of imaging.

Measure everything, 100x wider range

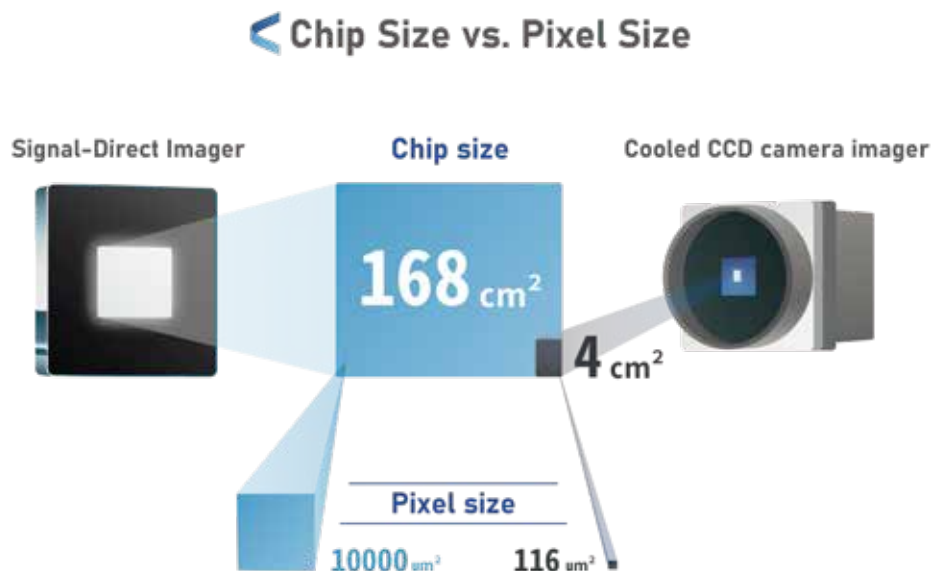

The sensor chip size impact overall signal collection. A larger chip can capture the signal from a wider area of the blot, which can improve the overall dynamic range and reduce the need for multiple captures or scanning. The Signal-Direct imager's large 168 cm^2 sensor captures the entire blot at once, avoiding the signal loss and artifacts that can occur with scanning smaller areas. In contrast, the Cooled CMOS camera's 4.0 cm^2 sensor is much smaller, which may require multiple images to be stitched together to capture a full blot.

The larger pixel area of the Signal-Direct Imager allows it to capture a significantly higher number of photons from the same signal source compared to a much smaller pixel on a standard CCD chip.

The Signal-Direct Imager pixel area ($100 \mu\text{m} \times 100 \mu\text{m}$) is approximately 277 to 1111 times larger than the Cooled CCD pixel area ($3 \text{ to } 6 \mu\text{m}$). This means it can collect hundreds to over a thousand times more photons, directly translating to a proportional increase in sensitivity.



The revolution in simplicity and power

Simplify your workflow.

From image capture to final results, our platform offers a seamless imaging experience. Perform instant data analysis then dig deeper with the robust Max-One software to get the answers you need, faster.

Leveraging an industry-leading to 6.0 OD dynamic range

It guarantees you capture true, publication-ready quantitative data by simultaneously detecting low-abundance targets and highly abundant controls without detector saturation.

This vast capacity provides a superior Signal-to-Noise Ratio for ECL while fully supporting the wider linear range of next-generation probes


Designed for Multi-User Labs Where Every Inch Matters.



The Signal-Direct Imager is engineered for efficiency, with a footprint so compact it's comparable to a standard laptop. That means more room for your other experiments, chemicals, and equipment.

Don't trade size for sensitivity: Get high-resolution imaging capability without dedicating a large room or extensive bench space.

Streamline Your Documentation: One Platform for All Your Results.

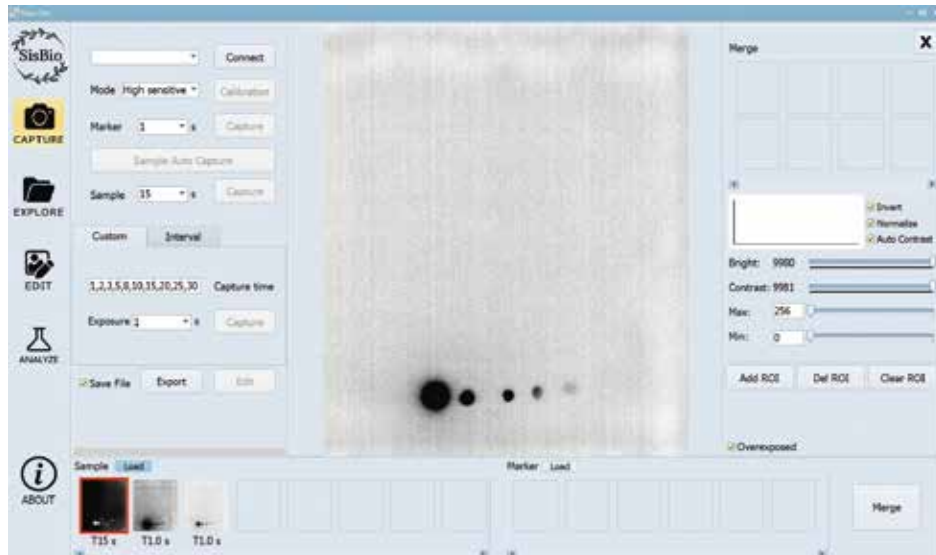

Essential Max-One Software Features

Dynamic Range Capture:

High Resolution and High Sensitivity in a Single Shot.


Why choose one mode when you can have both? Our software's dynamic capture technology automatically applies the optimal settings, ensuring you never miss a detail. It intelligently manages signal intensity to prevent saturation as much as possible, aiming to give you brilliant, high-quality images.

Optimal Contrast and Clarity for Every Shot.


Stop wasting time with multiple exposures. Our system automatically processes each image to bring out the maximum contrast and dynamic range possible, ensuring a perfect, data-rich image every time.

A faster workflow. A more beautiful image

The system provides a dedicated color protein marker capture capability for the most aesthetic and precise image documentation available

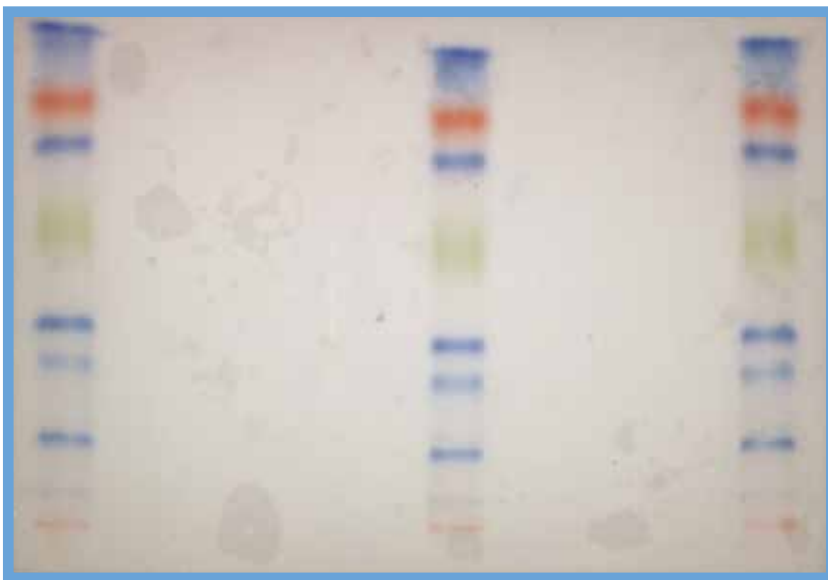

Achieve superior detection down to an exceptional 0.125 pg/ μ L (Dot Blot) in just 5 seconds using our Custom mode

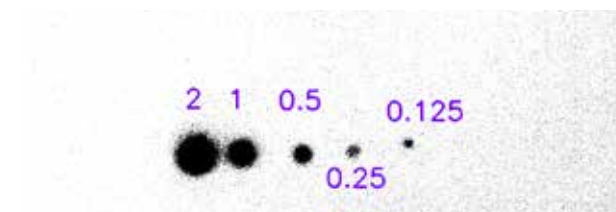

Mouse IgG HRP-conjugated Ab cat. no. : HAF007, bio-technie (R&D systems)
LumiFlash Femto cat. No. LF24-100, Visual protein
PVDF membrane (0.45um) E-BC-R266, Elabscience

Technical Specifications



Pixel Size	100 um x 100 um
CMOS sensor chip size	168 cm ²
Active Imaging Area (W x D)	13.7 cm x 11.7 cm (5.4 in. x 4.6 in.)
Full well capacity, electronics	3,150,000 (Maximum)
Dynamic range	6.0 OD
Operation Mode	Auto & Custom exposure
Data Format	TIFF/JPEG/BMP/PNG
Analysis software	"PC control, At least one available USB 3.0 (or newer/higher) port is required."
Dimension(W x D x H)	280 mm x 250 mm x 53mm(110 in x 98 in x 21 in)
Weight	Approximately 4 Kg (8.8 lb)
Power	100 ~ 240V, 50 ~60 Hz
Light Source	Chemiluminescence, epi LEDs (R/G/B)
Operation Temperature	15° C ~ 35° C
Operating Humidity	35% ~ 55% relative humidity (noncondensing)





SisBio Limited

Address: 4 F., No. 12, Ln. 345, Yangguang St., Neihu Dist.,
Taipei City 114713, Taiwan (R.O.C.)

Postal code 114713

Email: Sales@sisbio.com.tw

Service@sisbio.com.tw

Website: www.sisbio.com.tw

For Research Use Only. Not for use in diagnostic procedures

SisBio are trademarks of the company. All other trademarks
appearing in this brochure are the property of their respective owners